

■ SUCTION ASPIRATORS HIGH VACUUM, LOW AND HIGH FLOW


- 28222 TOBI - suction aspirator 230V - 50Hz
- 28224 SUPER TOBI - suction asp. 230V - 50Hz
- 28225 SUPER TOBI - suction asp. 110V - 60Hz

Portable suction aspirators, ideal for tracheotomy and small surgery. Vacuum continuously adjustable with vacuum indicator. Have unbreakable 1,000 ml standard bottle (2,000 ml optional) autoclavable at 120°C with safety float control valve to prevent overflow. Silicone connection tube. ABS plastic case. Made in Italy

Tobi: suction power 16 air litres/min
max suction -0.75 bar

Super tobi: suction power 40 air litres/min
max suction -0.80 bar

See spare accessories at page 82 and other accessories at page 84

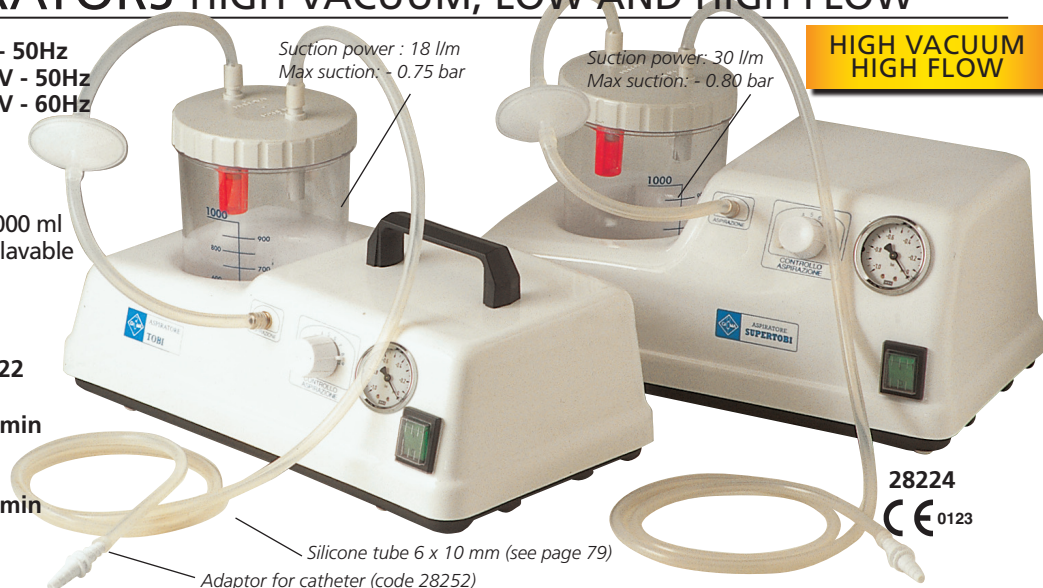

**HIGH VACUUM
HIGH FLOW**

28224
CE 0123

STANDARD ACCESSORIES

	Tobi	Super Tobi
Bottle 1000 ml with cover	1	1
Antibacterial Filter	1	1
Suction catheter	1	1
Silicon Tube set	1	1
User manual*	1	1

* User manual is provided in 5 languages: English, Spanish, French, German and Italian.

TECHNICAL SPECIFICATIONS

	28222 Tobi	28224/28225 Super Tobi
Operating voltage:	230 V - 50 Hz	230 V - 50 Hz or 110 V - 60 Hz
Power consumption:	184 W	106 W
Bottle capacity:	1000 ml	1000 ml
Flow (air litres/min):	18 l/min	40 l/min
Max suction:	-0.75 bar (563mm/Hg)	-0.80 bar (600mm/Hg)
Working time: minutes	20 ON / 40 OFF	120 ON / 60 OFF
Size (cm):	37x22x21 cm	37x22x21 cm
Weight:	3.5 kg	5.25 kg
Norms:	IEC 601-1	IEC 601-1

■ PROFESSIONAL ASPIRATORS WITH METAL BOX - HIGH VACUUM

- 28244 ASPEED - suction - 230 V - single pump
- 28245 ASPEED PROFESSIONAL - suction - 230 V - double pump
- 28246 ASPEED PROFESSIONAL - suction - 110 V - double pump

For professional surgical aspiration. It requires no maintenance or lubrication and is ideal for clinical use. The piston-type electric aspirator gives high performance and great durability. It is equipped with a protective thermal cut-out relay. It is also equipped with motor-protector cap that totally prevents aspirated liquids or secretions from reaching and damaging the vacuum pump.

- 28255 HYDROPHOBIC 99% BACTERIAL FILTER
- 28258 AUTOCLAVABLE JAR 1l with cover - since 2007
- 28259 DISPOSABLE SUCTION LINER
- 25850 SILICONE TUBE 6 X 12 mm - min. order 30 m
- 28254 ADAPTOR to connect catheter to silicone tube

STANDARD ACCESSORIES

Autoclavable polycarbonate jar 1000 cc with safety valve (overflow protection)
Disposable suction liner
99% Antibacterial hydrophobic filter
Sterile cannula
Sterile manual regulator
Set of atoxic sterilizable silicone tubes
Power Cable
User Manual (English, Italian, French, German, Spanish)



CE 0434
EN ISO 10079-1

TECHNICAL SPECIFICATIONS

	28244	28245/28246
Operating Voltage:	230V-50 Hz	230V-50 Hz/110V-60 Hz
Bottle capacity:	1 l	1 l
Flow:	15 air l/min	22 air l/min
Max suction		- 0.85 bar
Size:		18.5x33x21 cm
Weight:	3.5 kg	4.5 kg
Material:		Metal case
Sound level:	55 dBA	65 dBA
Norms:		Class I Type B



28258
+
28259