


VEGA SUCTION UNITS, HIGH VACUUM, LOW AND HIGH FLOW



**HIGH VACUUM
HIGH FLOW**

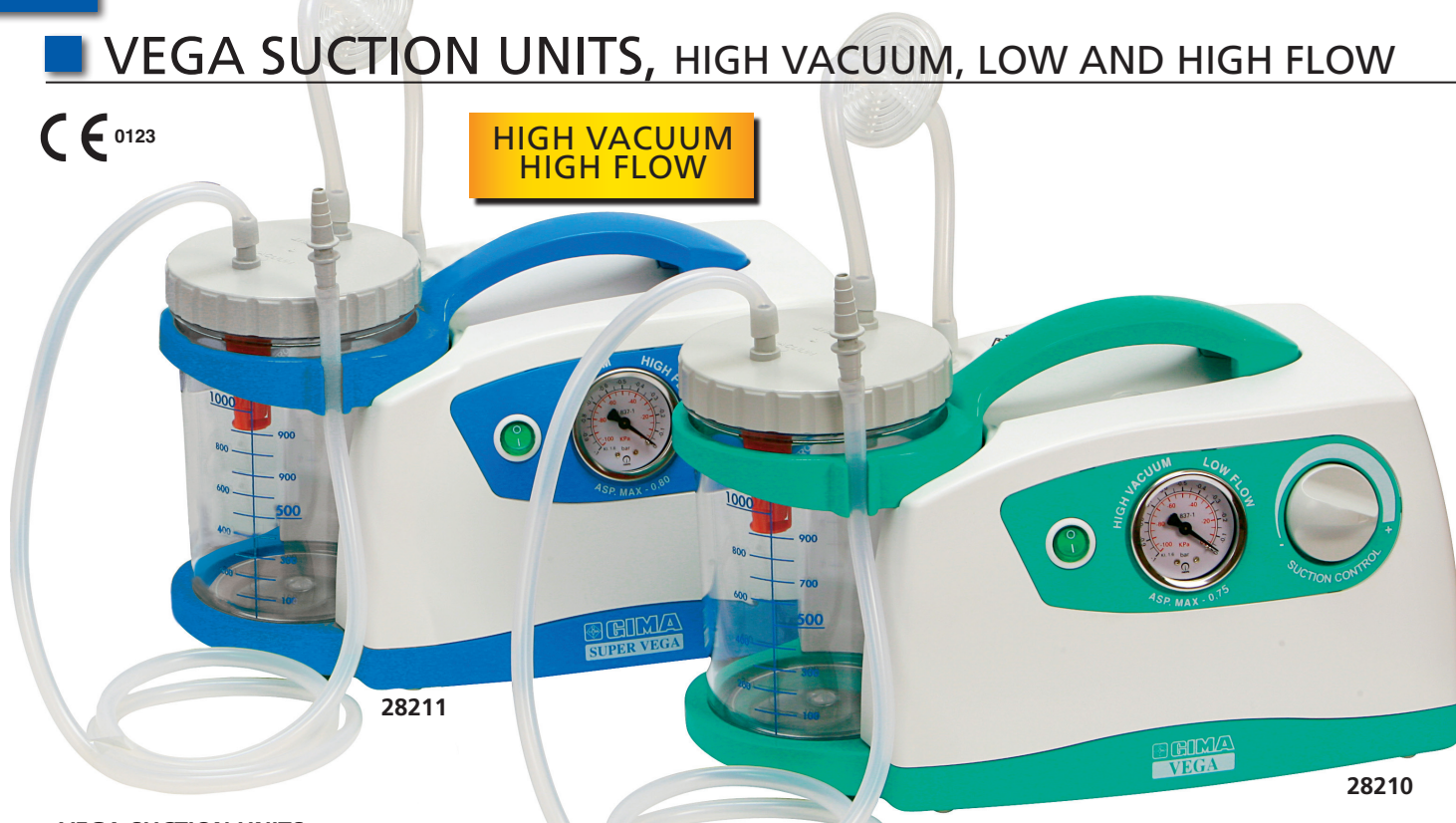

28211

28210

• VEGA SUCTION UNITS


A line of stylish, compact and reliable aspirators for small surgery. Silent, high flow and maintenance - free pumps allow long product life and durability.

Recommended for nasal aspiration, oral aspiration, tracheal aspiration of body liquids (mucus, catarrh or blood) in adults or children. Equipped with vacuum regulator continuously adjustable, vacuum indicator, unbreakable 1,000 ml or 2,000 ml bottle, autoclavable at 120°C, with safety float control valve to prevent overflow.

Manufactured in electric and thermal insulated plastic body.

Autoclavable and disposable collection system: In addition to standard autoclavable jar is available a disposable collection system composed by a rigid reusable container and a polyethylene disposable liner with cover, ermetically sealed, and hydrophobic, antireflux, bacterial filter operating as an overflow valve.

	28209 VEGA UNO	28210 VEGA	28211 SUPERVEGA
Operating Voltage:	230 V - 50 Hz 110 V - 60 Hz (on request)	230 V - 50 Hz 110 V - 60 Hz (on request)	230 V - 50 Hz 110 V - 60 Hz (on request)
Power consumption:	184 W	184 W	106 W
Bottle capacity:	1000 ml	1000 ml	1000 ml or 2000 ml
Flow (air litres/min):	15 l/min	16 l/min	40 l/min
Max suction:	-0.75 bar (563 mm/Hg)	-0.75 bar (563mm/Hg)	-0.80 bar (600mm/Hg)
Temporary use:	20/40 min on/off	30/30 min on/off	120/60 min on/off
Size:	25x19x16.5 cm	35x21x18 cm	35x21x18 cm
Weight:	2.2 kg	2.5 kg	3.6 kg
Norms:	93/42/ECC	93/42/ECC	93/42/ECC
Origin:	Italy	Italy	Italy

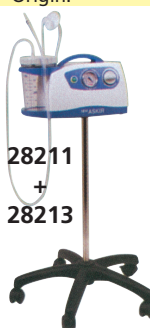

28211
+
28213

28252
Adaptor to connect
catheter to silicone
tube

GIMA code	VEGA SUCTION ASPIRATOR	Operating power	Flow
28209	Vega Uno	230 V - 50 Hz	15 l/m
28210	Vega	230 V - 50 Hz	16 l/m
28211	Super Vega 1l	230 V - 50 Hz	40 l/m
28212	Super Vega 2l	230 V - 50 Hz	40 l/m
-	Super Vega 1l - on request	110 V - 60 Hz	40 l/m



GIMA code	SPARE ACCESSORIES FOR VEGA AND TOBI LINE
28229	Hydrophobic, 99% antibacterial filter
28226	Autoclavable jar 1 l with cover and overflow valve system
28236	Autoclavable jar 2 l with cover and overflow valve system
28227	Cover and overflow valve system for 28226 and 28236
28273	Disposable liner 2 l (cover + bag)
28272	Container for disposable liner
25845	Silicone tube 6x10 mm - minimum order 30 m
28252	Adaptor to connect catheter to silicone tube
	Suction catheter - see at page
28213	Trolley for Vega and Super Vega

STANDARD ACCESSORIES

Autoclavable jar (1 or 2 l) with overflow valve system
Antibacterial and hydrophobic filter
Set of silicone tubes 6x10 mm
Conical connector for probes
5 languages user manual (English, French, German, Italian, Spanish)



1 l polycarbonate
autoclavable jar
with overflow
safety valve

28209